

## STOCK *Nr2e3<sup>rd7</sup>/J*

Stock No: 002139

Protocol 37734: Separated PCR Assay - Nr2e3<rd7> Alternate1

Version 1.0

### Notes

The genotyping protocol(s) presented here have been optimized for reagents and conditions used by The Jackson Laboratory (JAX). To genotype animals, JAX recommends researchers validate the assay independently upon receipt of animals into their facility. Reaction cycling temperature and times may require additional optimization based on the specific genotyping reagents used.

### Expected Results

Mutant = 149 bp

Heterozygote = 149 bp and 212 bp

Wild type = 212 bp

[chr9:59948296-59948507](#) 212bp TGATGTCACCAGCAATGACC CATGTTTCCCAAGCTGACT

### JAX Protocol

#### Protocol Primers

PRIMER	5' LABEL	SEQUENCE 5' → 3'	3' LABEL	PRIMER TYPE	REACTION	NOTE
51784		CAT GTT TCC CCA AGC TGA CT		Common	A, B	
52033		TGA TGT CAC CAG CAA TGA CC		Wild type Forward	A	Wt F
52034		CTG TTC AGA CTA ATT TCT GTG CTG		Mutant Forward	B	Mut F

#### Reaction A

COMPONENT	FINAL CONCENTRATION
ddH <sub>2</sub> O	
Kapa 2G HS buffer	1.30 X
MgCl <sub>2</sub>	2.60 mM
dNTPS-kapa	0.26 mM
51784	0.50 uM
52033	0.50 uM
Glycerol	6.50 %
Dye	1.00 X
Kapa 2G HS taq polym	0.03 U/ul
DNA	

#### Cycling

STEP	TEMP °C	TIME	NOTE
1	94.0	--	
2	94.0	--	
3	65.0	--	-0.5 C per cycle decrease
4	68.0	--	
5		--	repeat steps 2-4 for 10 cycles (Touchdown)
6	94.0	--	
7	60.0	--	
8	72.0	--	
9		--	repeat steps 6-8 for 28 cycles
10	72.0	--	
11	10.0	--	hold

JAX uses a very high speed Taq (~1000 bp/sec), use cycling times recommended for your reagents.

JAX uses a 'touchdown' cycling protocol and therefore has not calculated the optimal annealing temperature for each set of primers.

## Reaction B


COMPONENT	FINAL CONCENTRATION
ddH2O	
Kapa 2G HS buffer	1.30 X
MgCl2	2.60 mM
dNTPS-kapa	0.26 mM
51784	0.50 uM
52034	0.50 uM
Glycerol	6.50 %
Dye	1.00 X
Kapa 2G HS taq polym	0.03 U/ul
DNA	



## Cycling

STEP	TEMP °C	TIME	NOTE
1	94.0	--	
2	94.0	--	
3	65.0	--	-0.5 C per cycle decrease
4	68.0	--	
5		--	repeat steps 2-4 for 10 cycles (Touchdown)
6	94.0	--	
7	60.0	--	
8	72.0	--	
9		--	repeat steps 6-8 for 28 cycles
10	72.0	--	
11	10.0	--	hold

JAX uses a very high speed Taq (~1000 bp/sec), use cycling times recommended for your reagents.

JAX uses a 'touchdown' cycling protocol and therefore has not calculated the optimal annealing temperature for each set of primers.

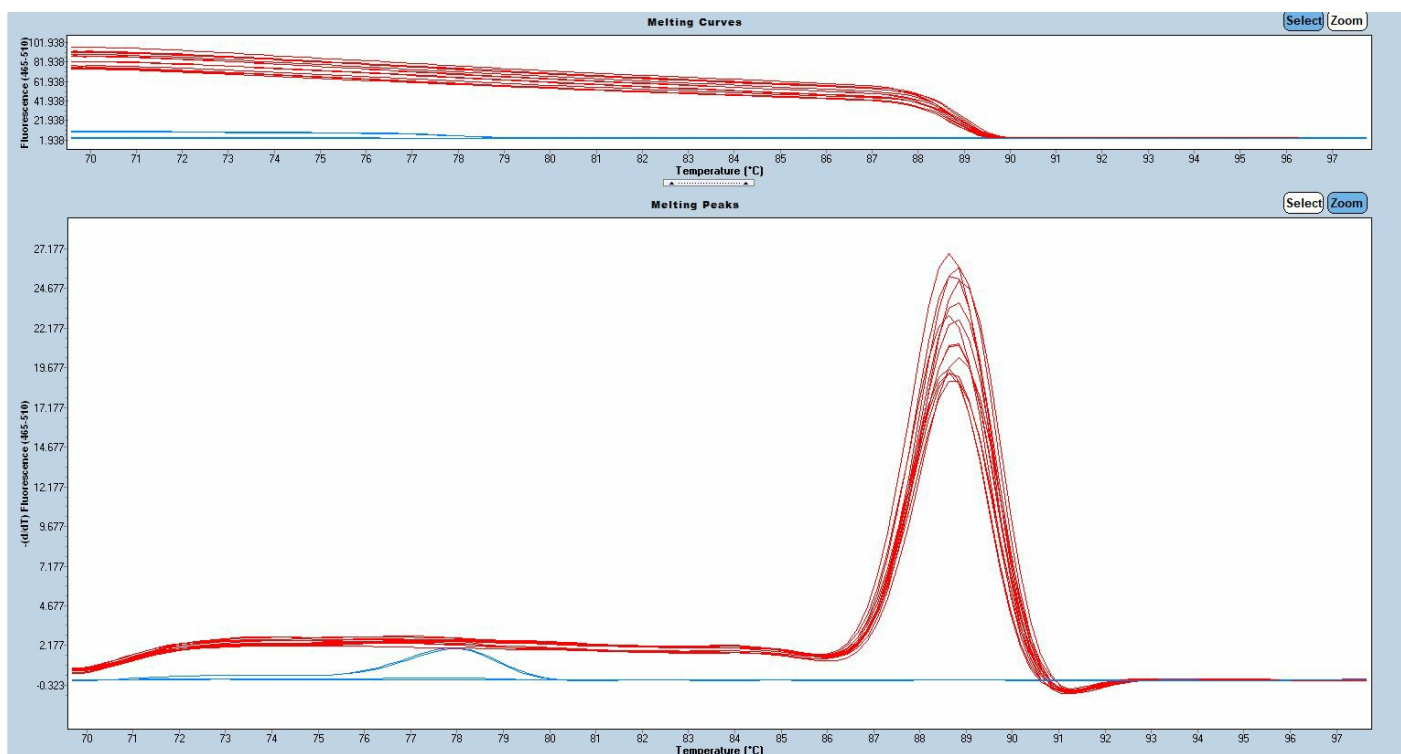

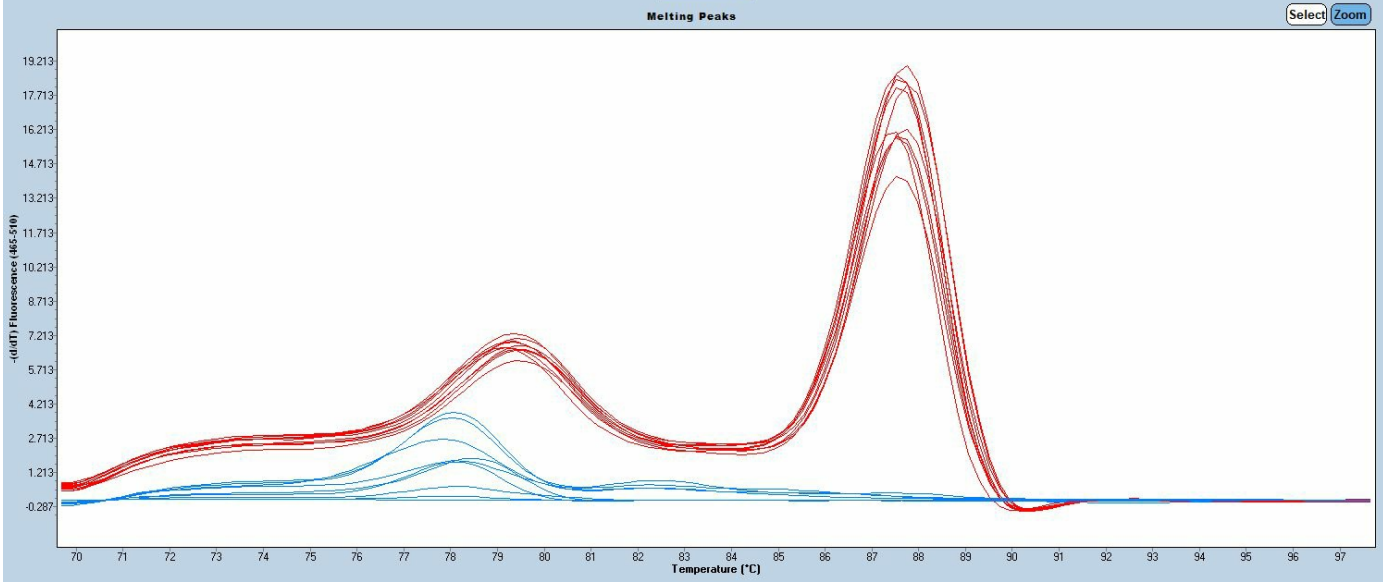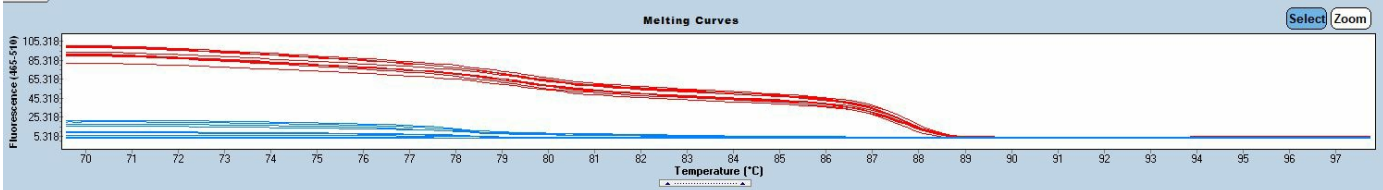

31255 Nr2e3rd7 Alt1

