

C57BL/6-Tg(Pmch-cre)1Rck/J

Stock No: 024034

Protocol 29382: Standard PCR Assay - Tg(Pmch-cre)1Rck alternate1

Version 1.2

Notes

This assay will NOT distinguish hemizygous from homozygous transgenic animals.

The genotyping protocol(s) presented here have been optimized for reagents and conditions used by The Jackson Laboratory (JAX). To genotype animals, JAX recommends researchers validate the assay independently upon receipt of animals into their facility. Reaction cycling temperature and times may require additional optimization based on the specific genotyping reagents used.

Expected Results

Transgene = ~250 bp

Internal positive control = 324 bp

JAX Protocol

Protocol Primers

PRIMER	5' LABEL	SEQUENCE 5' → 3'	3' LABEL	PRIMER TYPE	REACTION	NOTE
14314		CGG TTA TTC AAC TTG CAC CA		Transgene Reverse	A	
21173		CCG CAT CTT CTT GTG CAG T		Internal Positive Control Forward	A	
21175		ATC ACG TCC TCC ATC ATC C		Internal Positive Control Reverse	A	
22064		TCC AAT GCA CTC TTG TTT GG		Transgene Forward	A	

Reaction A


COMPONENT	FINAL CONCENTRATION
ddH ₂ O	
Kapa 2G HS buffer	1.30 X
MgCl ₂	2.60 mM
dNTP KAPA	0.26 mM
14314	0.50 uM
21173	0.50 uM
21175	0.50 uM
22064	0.50 uM
Glycerol	6.50 %
Dye	1.00 X
Kapa 2G HS taq polym	0.03 U/ul
DNA	

Cycling

STEP	TEMP °C	TIME	NOTE
1	94.0	--	
2	94.0	--	
3	65.0	--	-0.5 C per cycle decrease
4	68.0	--	
5		--	repeat steps 2-4 for 10 cycles (Touchdown)
6	94.0	--	
7	60.0	--	
8	72.0	--	
9		--	repeat steps 6-8 for 28 cycles
10	72.0	--	
11	10.0	--	hold

JAX uses a very high speed Taq (~1000 bp/sec), use cycling times recommended for your reagents.

JAX uses a 'touchdown' cycling protocol and therefore has not calculated the optimal annealing temperature for each set of primers.

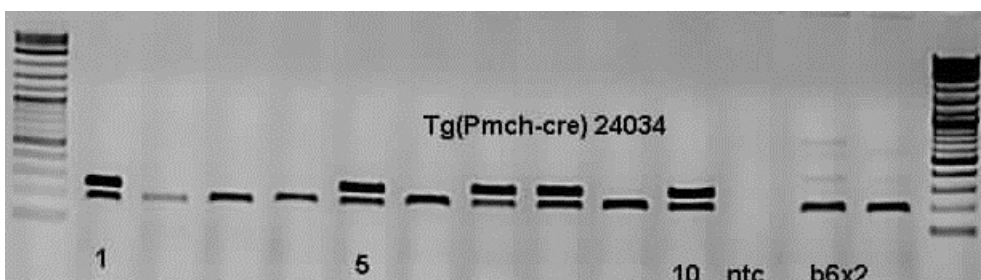

Melting Peaks

